无菌分装隔离系统 Sterile subpackage isolation system

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PRODUCT INTRODUCTION:

The filling aseptic isolator system consists of the following parts: aseptic isolation operation area, material transfer area, bag-in -bag-out pollutant collection and processor, pneumatic sealing system, control system, monitoring system, and a supporting VHP sterilization system.


According to the technological functions, it can be classified into: production type, weighing and sampling type, aseptic feeding type, aseptic filling type, etc. It is mainly used in the technological processes such as aseptic filling, weighing, sampling, filling, and packaging of pharmaceutical preparations, raw materials, and auxiliary materials.

Composition of the isolation cabin: receiving cabin, filling cabin, inner - cover - putting cabin, capping cabin, discharging cabin.

Main configurations: electric drum, vibrating filling machine, electronic scale, capping machine, VHP sterilization system, automatic barrel - holding machine, photoelectric - induction automatic sliding door, dust collector.

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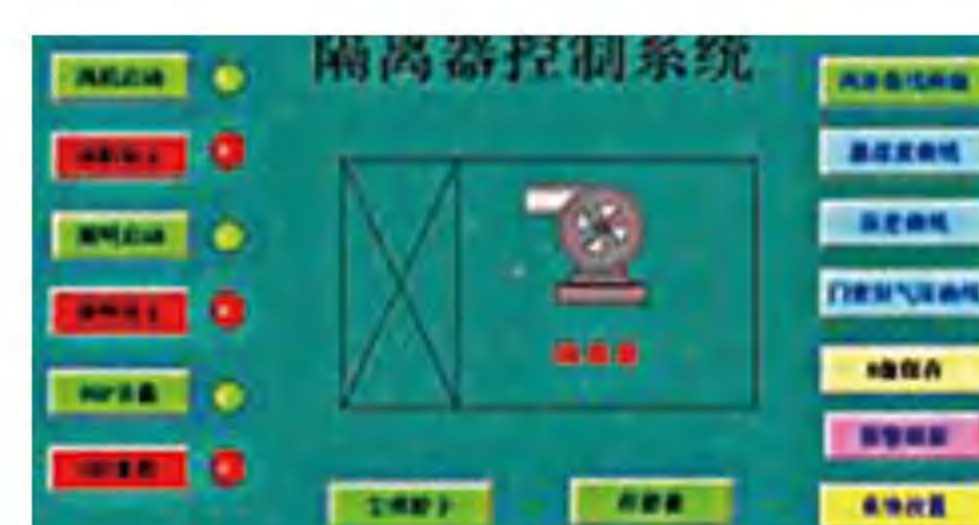TECHNICAL FEATURES:

Double - door transfer technology: Ensures that the clean space is not affected by material transfer and personnel operations, and solves the problem that the 100 - class cleanliness of traditional aseptic rooms cannot be continuously maintained.
Monitoring system: Online monitoring of operating parameters of the isolator, such as pressure, temperature, humidity, cleanliness, planktonic bacteria, wind speed, and exhaust volume, which meets the management requirements of GMP/GLP.
Control system: Automatically controlled by a microcomputer, with a large - screen touch - type human - machine interface, compliant with computerized validation software.
Adjustable positive and negative pressure: Can meet the protection requirements for toxic substances and prevent pollution.
Sealing method: Double - layer sealing to achieve zero leakage.
Sterilization modes: Two modes, ozone or vaporized hydrogen peroxide (VHPS), are available for configuration.

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CUSTOMER CASES: